



HAL
open science

Simple Fast Quantification of Cholecalciferol, 25-Hydroxyvitamin D and 1,25-Dihydroxyvitamin D in Adipose Tissue Using LC-HRMS/MS

Laurianne Bonnet, Marielle Margier, Ljubica Svilar, Charlène Couturier,
Emmanuelle Reboul, Jean-Charles Martin, Jean-Francois Landrier, Catherine
Defoort

► To cite this version:

Laurianne Bonnet, Marielle Margier, Ljubica Svilar, Charlène Couturier, Emmanuelle Reboul, et al.. Simple Fast Quantification of Cholecalciferol, 25-Hydroxyvitamin D and 1,25-Dihydroxyvitamin D in Adipose Tissue Using LC-HRMS/MS. *Nutrients*, 2019, 11 (9), pp.1977. 10.3390/nu11091977. hal-02487245

HAL Id: hal-02487245

<https://amu.hal.science/hal-02487245v1>

Submitted on 21 Feb 2020

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.


L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution 4.0 International License

Article

Simple Fast Quantification of Cholecalciferol, 25-Hydroxyvitamin D and 1,25-Dihydroxyvitamin D in Adipose Tissue Using LC-HRMS/MS

Laurianne Bonnet ^{1,†}, Marielle Margier ^{1,†}, Ljubica Svilar ^{2,†}, Charlene Couturier ¹, Emmanuelle Reboul ¹ , Jean-Charles Martin ^{1,2}, Jean-François Landrier ^{1,2} and Catherine Defoort ^{1,2,*}

¹ C2VN (Centre de Recherche CardioVasculaire et Nutrition), INRA, INSERM, Université d'Aix-Marseille, Faculté de Médecine, 27 Bd Jean Moulin, 13385 Marseille CEDEX 05, France

² CriBioM, Criblage Biologique Marseille, Faculté de médecine de la Timone, 13005 Marseille, France

* Correspondence: catherine.defoort@univ-amu.fr; Tel.: +33-4-91-32-42-82; Fax: +33-4-91-78-21-01

† These three authors contributed equally to this work.

Received: 17 June 2019; Accepted: 17 August 2019; Published: 22 August 2019



Abstract: Vitamin D metabolism is actively modulated in adipose tissue during obesity. To better investigate this process, we develop a specific LC-HRMS/MS method that can simultaneously quantify three vitamin D metabolites, i.e., cholecalciferol, 25-hydroxyvitamin D₃ (25(OH)D₃), and 1,25-dihydroxyvitamin D₃ (1,25(OH)₂D₃) in a complex matrix, such as mouse adipose tissue and plasma. The method uses pretreatment with liquid–liquid or solid–phase extraction followed by derivatization using Amplifex[®] reagents to improve metabolite stability and ionization efficiency. Here, the method is optimized by co-eluting stable isotope-labelled internal standards to calibrate each analogue and to spike biological samples. Intra-day and inter-day relative standard deviations were 0.8–6.0% and 2.0–14.4%, respectively for the three derivatized metabolites. The limits of quantification (LoQ) achieved with Amplifex[®] derivatization were 0.02 ng/mL, 0.19 ng/mL, and 0.78 ng/mL for 1,25(OH)₂D₃, 25(OH)D₃ and cholecalciferol, respectively. Now, for the first time, 1,25(OH)₂D₃ can be co-quantified with cholecalciferol and 25(OH)D₃ in mouse adipose tissue. This validated method is successfully applied to study the impact of obesity on vitamin D status in mice.

Keywords: vitamin D metabolite quantification; LC-MS/MS; adipose tissue; obesity

1. Introduction

Vitamin D (VD), or cholecalciferol, is a prohormone involved in calcium and phosphate homeostasis. VD insufficiency, or deficiency, is directly linked to well-known bone diseases including rickets, osteoporosis, and osteomalacia [1], and may be involved in a number of other pathophysiological disorders, including cardiovascular disease or type 2 diabetes. After endogenous synthesis or intestinal absorption [2], cholecalciferol is metabolized in the liver by 25-hydroxylases that produce 25-hydroxyvitamin D [3]. The major form of VD found in plasma is 25(OH)D₃, which is considered the best biomarker of VD status [4]. 25(OH)D can be hydroxylated in the kidney or in target tissues [5] by CYP27B1 into 1 α ,25-dihydroxyvitamin D, which is the active form of VD.

Several studies and a recent review [6,7] highlighted that adipose tissue is one of the main storage sites for cholecalciferol and 25(OH)D₃, and the mechanism of cholecalciferol and 25(OH)D₃ uptake by adipocytes was recently identified [8]. Obesity is commonly associated with vitamin D deficiency [9,10], but there is no clear picture or even consensus on the mechanism(s) underpinning vitamin D status in obesity. One often-cited mechanism is deposition and sequestration of vitamin D metabolites in

adipose tissue stores, but there is little experimental evidence to support this hypothesis. We recently demonstrated that cholecalciferol and 25(OH)D₃ storage could result from an active mechanism involving modulated expression of several cytochromes [11].

One of the most informative ways to explore VD-related health and metabolism is to simultaneously measure VD metabolites in key organs. Consequently, the frequency of routine analytical measurement of VD status has increased greatly [12], necessitating the development of reliable and accurate analytical methods to measure concentration levels of VD metabolites not only in plasma but also in relevant body tissues [13]. Numerous liquid chromatography-tandem mass spectrometry (LC-MS/MS) methods have been described for measuring VD metabolites individually or in combination [14–21], but few papers reported the absolute quantification of cholecalciferol and 25(OH)D₃ in adipose tissue by high-performance liquid chromatography tandem-mass spectrometry (HPLC-MS/MS) [22–24], and none of them include the detection of 1,25(OH)₂D₃ with satisfactory sensitivity [25,26].

Here, to address this gap, we set out to develop and validate a liquid-chromatography high-resolution tandem mass spectrometry (LC-HRMS/MS) method for simultaneous and accurate determination of cholecalciferol, 25(OH)D₃ and 1,25(OH)₂D₃ in mouse adipose tissue. This effort carried a number of challenges, chiefly (i) the low in-sample concentrations of these metabolites, especially 1,25(OH)₂D₃ which is biologically active at picomolar concentrations, and (ii) the small amounts of tissue samples available, especially samples from mouse studies. The need to increase ionization efficiency and sensitivity implied a sample pretreatment, so we opt to use Amplifex[®] as a derivatizing reagent applicable to diene analytes. We apply our validated method to quantify the three vitamin D metabolites in adipose tissue and plasma from C57BL/6J mice fed for 11 weeks on either a control diet (control, 10% energy from fat) or a high-fat diet (HF, 60% energy from fat) formulated to provide equivalent vitamin D intake in both groups.

2. Material and Methods

2.1. Materials

Cholecalciferol, 25(OH)D₃, 1,25(OH)₂D₃, dibasic potassium phosphate (K₂HPO₄), methyl tert-butyl ether (MTBE), and phosphate buffer saline (PBS) were purchased from Sigma–Aldrich. Amplifex[®] Diene reagent was purchased from Sciex (Chemistry and Consumables R&D, Framingham, MA). Isotopically-labelled d₃-cholecalciferol, d₃-25(OH)D₃, and d₃-1,25(OH)₂D₃ were purchased from Cambridge Isotope Laboratories (MA). These three labelled molecules were deuterated in position 6,19,19 and have been used as internal standards for the precision of the method. All vitamin D metabolites were stored at –80 °C in the dark. LC-MS-grade acetonitrile (AcN) and formic acid and HPLC-grade methanol (MeOH) and ethyl acetate (EtOAc) were obtained from Carlo Erba Reagents (Peypin, France). The water used was generated through a Direct-Q Ultrapure Water System from Millipore (Bedford, MA) with a specific resistance of 18.2 MΩ cm. Oasis HLB cartridges used for solid-phase extraction (3 CC/60 mg) were supplied by Waters (Guyancourt, France). Mouse samples were taken on euthanized animals (C57Bl/6J male mice) according to the local animal ethics committee protocols (APAFIS#2595-2016091911217758).

2.2. Instrumentation

Chromatographic separation was carried out on a Dionex UltiMate 3000 system (Thermo Fisher Scientific, Waltham, MA) configured with a rapid separation (RS) pump (LPG-3400 RS), an autosampler (WPS-3000 TRS) and a column compartment (TCC-3000 RS), all operated by Chromeleon 6.8 software. Samples were injected on a Hypersil GOLD C18 column (2.1 × 100 mm; Thermo Scientific, Les Ulis, France). Parallel reaction monitoring (PRM) and accurate mass measurements were performed on a Q-Exactive Plus mass spectrometer (Thermo Fisher Scientific, Bremen, Germany) with a heated electrospray ionization (H-ESI II) probe. Thermo Xcalibur 3.0.63 software was used for both instrument

setup and LC-MS system control during data acquisition, and for data treatment. Q Exactive Plus Tune 2.5 software was used for direct control of the mass spectrometer.

2.3. Chromatography Conditions

The autosampler tray was kept at 4 °C during the experiment, and column oven temperature was regulated at 40 °C. 5 µL of sample was injected onto the column. Mobile phase A consisted of 0.1% formic acid in water, and mobile phase B was 0.1% formic acid in AcN. Solvent flow rate was set to 400 µL/min. Starting mobile phase consisted of 30% of solvent B held for 4 min. A linear gradient was applied with B % increased to 65% for 6 min, held isocratic for 2 min, then increased to 100% of B for 4 min. Start conditions were reinstated in 2 min, making the total run time 18 min.

2.4. Mass Spectrometry Conditions

LC-HRMS/MS analyses were performed with external calibration in positive mode providing a mass precision lower than 5 ppm. The H-ESI probe was kept at 310 °C. Transfer capillary temperature was kept at 320 °C. Spray voltage was set at 3500 V and the S-lens RF level was set at 55 V. Sheath and auxiliary gas were held at 30 and 8 (arbitrary units).

PRM was set up in MS/MS acquisition mode, and full higher-energy collision dissociation (HCD) spectra were acquired. For each compound, the most abundant ion was selected as precursor ion and isolated in a 2 U_ma window in a specified time segment and fragmented under previously optimized collision energies (Table 1). Resolving power was set to 35,000 full width at half maximum (FWHM) at *m/z* 200 and the Auto Gain Control (AGC) was set to 2×10^5 . Injection time was set to a maximum of 100 ms in order to optimize the analytical cycle time.

Table 1. Precursor and product ions *m/z* values, retention times, and normalized collision energies for Amplifex-derivatized cholecalciferol, 25(OH)D₃, 1,25(OH)₂D₃, and the corresponding isotopically-labeled molecule for absolute quantification.

Compound Name	Precursor Ion	Precursor <i>m/z</i>	Product <i>m/z</i>	Retention Time (min)	Collision Energy (%)
Cholecalciferol-Amplifex	[M] ⁺	716.50	657.43	13.2 + 13.5	21
25(OH)D ₃ -Amplifex	[M] ⁺	732.51	673.43	9.7	21
1,25(OH) ₂ D ₃ -Amplifex	[M] ⁺	748.51	689.42	8.6	21
d ₃ -cholecalciferol-Amplifex	[M] ⁺	719.50	660.45	13.2 + 13.5	21
d ₃ -25(OH)D ₃ -Amplifex	[M] ⁺	735.50	676.45	9.7	21
d ₃ -1,25(OH) ₂ D ₃ -Amplifex	[M] ⁺	751.50	692.44	8.6	21

2.5. Preparation of stock and working solutions of analytical and deuterated standards (IS)

A working solution of deuterated analytes was prepared at 0.02 ng/mL of each internal standard (IS) (d₃-cholecalciferol, d₃-25(OH)D₃ and d₃-1,25(OH)₂D₃) in ethanol. Stock solutions of cholecalciferol, 25(OH)D₃ and 1,25(OH)₂D₃ standards were prepared at 100, 50 and 10 ng/mL, respectively in ethanol and stored at −80 °C in the dark. Calibration curves were prepared by serial dilution of the three analyte stock solutions to obtain calibration standards from 0 to 50 ng/mL. After 1.5 µL of the working solution of deuterated analytes was added to each dilution, derivatization was performed using Amplifex[®]. A batch of Quality Control (QC) samples was produced for method validation.

2.6. Animals, Diets, and Experiments

Six-week-old male C57BL/6J mice were obtained from Janvier Labs (Le Genest-Saint-Isle, France) and fed ad libitum with standard chow (maintenance diet A04, Safe Diets, Augy, France) during a one week acclimatization period with ad libitum access to drinking water, housed at 22 °C under a 12 h/12 h light/dark cycle at 20% relative humidity. The mice were then divided into two groups: a control-diet group (control: 10% energy from lipids, *n* = 10) and a high-fat-diet group (HF: 60% energy

from lipids, $n = 10$) (TestDiet, London, UK). Blood and adipose tissue were collected after 11 weeks on either the control or high-fat diet.

Composition of the experimental diet is detailed in Supplementary Table S1. Morphology (body weight and adiposity index) and food intake (energy intake and vitamin D intake) parameters were also controlled and measured (Supplementary Tables S2 and S3, respectively). Further details can be found in Bonnet et al. (2019) [11].

2.7. Plasma Sample Preparation

At the end of the 11-week diet period, the mice were fasted overnight and blood was collected by cardiac puncture under anesthesia. Plasma was obtained after centrifugation at $3000\times g$ for 15 min at $4\text{ }^{\circ}\text{C}$, and stored at $-80\text{ }^{\circ}\text{C}$. The plasma preparation procedure was inspired by Wang et al. [27]. As vitamin D and its metabolites are light-sensitive, the extraction procedure was conducted under low light. After thawing on ice, mouse plasma was centrifuged at $13,600\times g$ for 15 min at $4\text{ }^{\circ}\text{C}$, then $100\text{ }\mu\text{L}$ of each sample was transferred to a glass tube with $10\text{ }\mu\text{L}$ of deuterated working solution. Proteins were precipitated by adding $500\text{ }\mu\text{L}$ of AcN, vortex-mixed, and centrifuged at $3000\times g$ for 10 min at $4\text{ }^{\circ}\text{C}$. The supernatant was transferred to another glass tube, then the volume was reduced to one half under a nitrogen stream, and $500\text{ }\mu\text{L}$ of EtOAc was added to the solution for liquid-liquid extraction. After shaking vigorously for 10 min, samples were centrifuged at $590\times g$ for 20 min at $4\text{ }^{\circ}\text{C}$, and the upper organic layer was transferred to a glass tube and reduced under nitrogen stream. The plasma samples were then derivatized with Amplifex[®] Diene reagent.

For validation, a set of plasmas were pooled together and treated as described above.

2.8. Adipose Tissue Preparation

The animals were sacrificed by cervical dislocation under anesthesia, and the epididymal white adipose tissue (eWAT) was collected, weighed, snap-frozen in liquid nitrogen, and stored at $-80\text{ }^{\circ}\text{C}$. The adipose tissue preparation procedure was inspired by Lipkie et al. [24]. After thawing on ice, 50 mg of tissue was ground with 1 mL of PBS at 30 mvt/s or 15 min (Retsch, Eragny-sur-Oise, France). Tissue homogenate was transferred into a glass tube and spiked with $25\text{ }\mu\text{L}$ of deuterated stock solution. AcN (1 mL) was added, and the mixture was vortexed for 5 min then centrifuged at $6000\times g$ for 5 min at $4\text{ }^{\circ}\text{C}$. After addition of 1 mL of MTBE, the mixture was vortexed for 5 min, then centrifuged at $6000\times g$ for 5 min at $4\text{ }^{\circ}\text{C}$, and the upper organic layer was collected in another tube. The same extraction was repeated twice more with additions of $300\text{ }\mu\text{L}$ of PBS, 1 mL of AcN (vortex + centrifugation) and 1 mL of MTBE (vortex + centrifugation). The supernatants were pooled then dried under nitrogen. Oasis HLB cartridges were prepared with 3 mL each of EtOAc, MeOH, and H_2O . The samples were reconstituted with 1 mL of MeOH and 1 mL of K_2HPO_4 (0.4 M) and loaded on the cartridge. The cartridge was washed with 3 mL of H_2O and 2 mL of 70% MeOH before being dried for 2 min under vacuum. Tips were washed with AcN and analytes were eluted with 1.5 mL of AcN and dried under nitrogen. Samples were then derivatized with Amplifex[®] reagent.

For validation, a set of adipose tissues were pooled together and treated as described above.

2.9. Amplifex Derivatization

A one-step derivatization using Amplifex[®] Diene as reagent was employed to improve the ionization efficiency of the metabolites [17]. The dried sample was added with $30\text{ }\mu\text{L}$ of Amplifex[®] solution, vortexed for 15 s, then incubated for 30 min at room temperature. Next, $30\text{ }\mu\text{L}$ of deionized water was added to the sample, vortexed for 15 s, and transferred into vials for LC injection.

2.10. Method Validation

Our method was validated according to FDA criteria adapted from Shah et al. [28], and we evaluated specificity, selectivity, linearity, sensitivity, accuracy, precision, stability, and dilution tests.

2.10.1. Linearity and Limits of Quantification (LoQ)

In order to assess the linearity of the quantification methods, we analyzed increasing amounts of analytical-standard solution spiked with the deuterated pool of the three corresponding analytes. For each derivatization reagent and each analyte, i.e., cholecalciferol, 25(OH)D₃ and 1,25(OH)₂D₃, calibration curves were plotted with peak area ratio of vitamin D metabolite to the respective internal standard versus a range of analyte concentrations. For cholecalciferol, these two peaks represent diastereomers 13.2- and 13.5-min retention times. As the tandem-mass spectra showed the same pattern (same product ion and relative intensities), we summed their peak areas for the quantification.

Five concentrations of standards were chosen to cover the expected concentration range in the study. The series of concentrations were 0.20–50 ng/mL for cholecalciferol, 0.10–12.5 ng/mL for 25(OH)D₃, and 0.020–2.5 ng/mL for 1,25(OH)₂D₃. We processed five replicates for each concentration level.

The LoQ was set as the concentration at which precision was less than 20% of the relative standard deviation, with an accuracy of between 80% and 120% of the theoretical value [29].

2.10.2. Accuracy and Precision

Replicated samples at three concentration levels were analyzed in separate runs to determine intra-day and inter-day precision and accuracy.

Intra-day accuracy and precision were calculated by processing five replicates at three concentration levels. Inter-day accuracy and precision were determined by analyzing five replicate samples in 3 batches ($n = 15$) at three concentration levels.

Precision of the assay was estimated by the coefficient of variation for each concentration level. Accuracy was represented as the relative error to the nominal concentrations (RE %), calculated as ratio of the deviation to the theoretical value $\times 100$.

3. Statistical Analysis

The data are reported as means \pm SEM. Significant differences were determined using a Student's *t*-test, with the threshold for statistical significance set at $p < 0.05$.

4. Results and Discussion

4.1. Method Validation

4.1.1. LC–HRMS/MS

The optimal separation was achieved with linear gradient conditions using AcN/H₂O (with 0.1% FA) as mobile phase on a C18 column, as shown in Figure 1.

The three vitamin D-Amplifex[®] derivatives were chromatographically resolved from each other (Peak Resolution: $R_s > 1.5$). Tandem mass spectrometry analysis combined with satisfactory chromatographic separation ensured high specificity and selectivity of the method (Supplementary Figures S1 and S2).

The recently-developed reagent Amplifex is a 4-phenyl-1,2,4-triazoline-3,5-dione (PTAD) base with a naturally ionized quaternary ammonium group. It further enhances ionization efficiency compared to PTAD and thus achieves lower detection limits. This may explain the higher sensitivity of the Multiple Reaction Monitoring (MRM) methods for quantification of Amplifex-derivatized VD compared to PTAD-derivatized VD. As expected, the Amplifex reagent gave a lower LoQ than PTAD and thus emerges as the preferred method for derivatization. The scope of LC-MS/MS methods published to date has been limited to the quantification of one or two vitamin D metabolites in a single run. There is now a consensus that derivatization techniques can maximize both selectivity and sensitivity, and here we show the value of this derivatization method for enabling simultaneous and accurate LC–HRMS/MS quantification of cholecalciferol, 25(OH)D₃, and 1,25(OH)₂D₃ in plasma and adipose tissue [16,17,19,20].

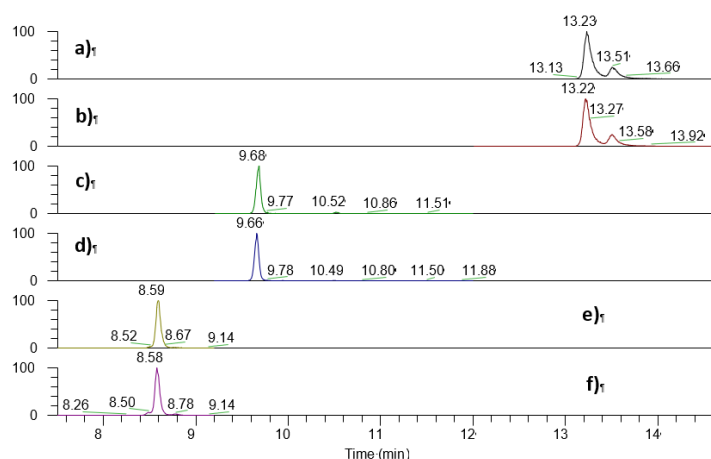


Figure 1. Extracted Ion Chromatography of standards mixed solution derivatized with Amplifex reagent. (a) m/z 657.43 extracted from the MS/MS spectrum of the m/z 716.50 ion of cholecalciferol, (b) m/z 660.45 extracted from the MS/MS spectrum of the m/z 719.50 ion of d3-cholecalciferol, (c) m/z 673.43 extracted from MS/MS spectrum of the m/z 732.51 ion of 25(OH)D3, (d) m/z 676.45 extracted from MS/MS spectrum of the m/z 735.50 ion of d3-25(OH)D3, (e) m/z 689.42 extracted from MS/MS spectrum of the m/z 748.51 ion of 1,25(OH)2D3 and (f) m/z 692.44 extracted from MS/MS spectrum of the m/z 751.50 ion of d3-1,25(OH)2D3.

4.1.2. Linearity, LoQ and Calibration Curves

The data were fitted using linear regression. Table 2 reports linear range and LoQ for the derivatized metabolites. Regression analysis showed good linearity with correlation coefficients (r^2) >0.993 for the three analytes.

Table 2. Linear range, limits of quantification (LoQ) and RE % for the 3 metabolites derivatized with Amplifex[®]. RE was determined for LoQ concentration (RE %: Relative error to the nominal concentrations).

Analytes	Linear Range ng/mL	LoQ ng/mL	RE %
Cholecalciferol	0.20–50	0.78	−5.3
25(OH)D3	0.10–12.5	0.19	4.0
1,25(OH)2D3	0.02–2.5	0.02	17.3

The range of linearity range was considered as the working range.

4.1.3. Precision and Accuracy

Table 3 summarizes the precision (%CV) and accuracy (%RE) of the method according to derivatization procedure for simultaneous quantification of all three vitamin D metabolites. Using Amplifex[®] derivatization, the intra-assay variation ranged between 0.8% and 6% depending on the vitamin D metabolite and concentration considered, and the inter-assay variation was between 2% and 14.4%. The RE of the three tested concentrations in three independent assays were within $\pm 3\%$, $\pm 9\%$ and $\pm 14\%$ for cholecalciferol, 25(OH)D3 and 1,25(OH)2D3, respectively. These results indicate good reproducibility and accuracy of the method. The %CV and %RE are all <15%.

4.1.4. Dilution Test and Stability

The dilution test with high-concentration overloaded plasma samples showed between 83% and 110% recovery (data not shown).

We tested the stability of derivatized vitamin D metabolites. The Amplifex-derivatized metabolites showed a signal decrease after 24 h in the autosampler (5 °C) from $-15.3 \pm 0.5\%$, $-19.7 \pm 2.0\%$, $-42.6 \pm 2.0\%$ for cholecalciferol, 25(OH)D3 and 1,25(OH)2D3, respectively (Data not shown). These stability results fit with what we expected based on Amplifex[®] Diene reagent guidance for users.

Table 3. Precision and accuracy of cholecalciferol, 25(OH)D3 and 1,25(OH)2D3 measurements.

Compound and Concentration	Intra-Assay (<i>n</i> = 5) ^a			Inter-assay (<i>n</i> = 3) ^b		
	Measured ^c	%RE ^d	%CV ^e	Measured ^c	%RE ^d	%CV ^e
Cholecalciferol						
1.56 ng/mL	1.53 ± 0.01	−1.7	1.0	1.57 ± 0.09	0.9	12.7
3.13 ng/mL	3.31 ± 0.06	6.2	2.6	3.09 ± 0.28	−1.2	9.0
12.5 ng/mL	12.9 ± 0.13	3.2	1.3	12.1 ± 0.25	−2.9	2.0
25(OH)D3						
0.78 ng/mL	0.81 ± 0.03	4.1	6.0	0.74 ± 0.09	9.3	11.1
1.56 ng/mL	1.70 ± 0.04	9.0	3.7	1.49 ± 0.10	−4.2	7.1
6.25 ng/mL	6.36 ± 0.12	1.8	2.7	5.83 ± 0.18	−6.8	3.1
1,25(OH)2D3						
0.16 ng/mL	0.15 ± 0.001	−6.3	1.5	0.18 ± 0.01	14.4	13.0
0.31 ng/mL	0.32 ± 0.002	5.1	0.8	0.31 ± 0.04	0.1	14.4
1.25 ng/mL	1.32 ± 0.02	5.8	1.9	1.23 ± 0.03	−1.6	2.8

^a intra-assay precision was calculated for 5 replicates measured in a single run; ^b inter-assay precision was calculated from 3 independent assays. ^c Means ± SD of measured concentrations. ^d Percent of relative error (%RE) is defined as the ratio of the deviation to the theoretical value ×100. ^e Coefficient of variance is defined as the ratio of the standard deviation to the mean value ×100.

4.2. Impact of High-Fat Diet on Vitamin D Metabolites in Plasma and in Adipose Tissue

After validation of the method, we quantified the three vitamin D metabolites in plasma and adipose tissue from the mouse study (Supplementary Figure S3).

In plasma, after 11 weeks of HF diet, compared to control diet, cholecalciferol concentration significantly decreased whereas 25(OH)D3 concentration increased and 1,25(OH)2D3 concentration remained unchanged. In eWAT, HF diet decreased cholecalciferol and 1,25(OH)2D3 concentrations but had no effect on 25(OH)D3 concentrations (Table 4). However, 25(OH)D3 was stronger in quantity terms (3.2-fold increase). These results are in agreement with Carrelli et al. [30], supporting the hypothesis that the large amount of adipose tissue in obese individuals serves as a reservoir for vitamin D. Here, our findings push the hypothesis further, with the quantification of 1,25(OH)2D3 concentration pointing to an active role of adipose tissue in the modulation of vitamin D metabolism during obesity. Note too that 1,25(OH)2D3 was also higher in quantity in obese mice compared to control mice (Table 4).

Table 4. Concentrations and quantities of the three Amplifex-derivatized vitamin D₃ metabolites in mouse plasma and adipose tissue after 11 weeks of diet. The quantities in adipose tissue were calculated as concentration × mass of epididymal white adipose tissue (eWAT).

	Cholecalciferol	25(OH)D3	1,25(OH)2D3
Plasma concentration			
Control	2.6 ± 0.1 ng/mL	38.4 ± 0.9 ng/mL	0.24 ± 0.02 ng/mL
HF	1.3 ± 0.1 ng/mL *	42.7 ± 1.4 ng/mL *	0.24 ± 0.01 ng/mL
Adipose tissue concentration			
Control	332 ± 25 ng/g	51 ± 4 ng/g	45 ± 6 ng/g
HF	236 ± 31 ng/g *	62 ± 8 ng/g	26 ± 6 ng/g *
Adipose tissue quantity			
Control	284 ± 29 ng	44 ± 5 ng	40 ± 6 ng
HF	548 ± 61 ng *	140 ± 12 ng *	61 ± 14 ng *

Control diet: *n* = 10, HF diet, high-fat diet: *n* = 10. Values are presented as means ± SEM. Significant differences were determined by a Student's *t*-test * *p* < 0.05.

5. Conclusions

The use of a recently developed derivatization reagent Amplifex[®] affords dramatically enhanced ionization efficiency, thus achieving lower detection limits and making it possible to detect picomolar-concentration metabolites. We have developed and validated a powerful LC-HRMS/MS technique to simultaneously quantify key analytes of interest in complex samples other than plasma. Here we used the method to quantify 1,25(OH)₂D₃ along with cholecalciferol and 25(OH)D₃ in adipose tissue to extend research into the hypothesis of vitamin D metabolite storage and production, especially in obesity settings [11].

Supplementary Materials: The following are available online at <http://www.mdpi.com/2072-6643/11/9/1977/s1>, Figure S1: MS/MS spectrum of 1,25(OH)₂D₃, 25(OH)D₃ and cholecalciferol derivatized with Amplifex reagent and their respective deuterated forms, Figure S2: Derivatization reaction, Figure S3: Extrated Ion Chromatography of plasma (A) and adipose tissue (B) sample derivatized with Amplifex reagent (a) *m/z* 657.43 extracted from the MS/MS spectrum of *m/z* 716.50 ion of cholecalciferol-Amplifex, (b) *m/z* 660.45 extracted from the MS/MS spectrum of *m/z* 719.50 ion of d₃-cholecalciferol-Amplifex, (c) *m/z* 673.43 extracted from MS/MS spectrum of *m/z* 732.51 ion of 25(OH)D₃-Amplifex, (d) *m/z* 676.45 extracted from MS/MS spectrum of *m/z* 735.50 ion of d₃-25(OH)D₃-Amplifex, (e) *m/z* 689.43 extracted from MS/MS spectrum of *m/z* 748.51 ion of 1,25(OH)₂D₃-Amplifex and (f) *m/z* 692.44 extracted from MS/MS spectrum of *m/z* 751.50 ion of d₃-1,25(OH)₂D₃-Amplifex, Table S1: Experimental diet compositions, Table S2: Mice morphologic parameters, Table S3: Food intake parameters.

Author Contributions: C.D., L.S., J.-F.L. conceived and designed the experiments; L.B. and M.M. performed the experiments; L.B., M.M., L.S. and C.D. analyzed the data; C.C., J.-C.M. contributed to reagents/materials/analysis tools; C.D. and L.B. wrote the paper with the contribution of E.R. and J.-F.L.

Funding: This research received no external funding.

Conflicts of Interest: The authors declare no conflicts of interest.

References

1. Plum, L.A.; DeLuca, H.F.; Vitamin, D. Disease and therapeutic opportunities. *Nat. Rev. Drug Discov.* **2010**, *9*, 941–955. [[CrossRef](#)] [[PubMed](#)]
2. Reboul, E.; Goncalves, A.; Comera, C.; Bott, R.; Nowicki, M.; Landrier, J.F.; Jourdheuil-Rahmani, D.; Dufour, C.; Collet, X.; Borel, P.; et al. Intestinal absorption is not a simple passive diffusion: evidences for involvement of cholesterol transporters. *Mol. Nutr. Food Res.* **2011**, *55*, 691–702. [[CrossRef](#)] [[PubMed](#)]
3. Schuster, I. Cytochromes P450 are essential players in the vitamin D signaling system. *Biochim. Biophys. Acta* **2011**, *1814*, 186–199. [[CrossRef](#)] [[PubMed](#)]
4. Holick, M.F.; Vitamin, D. Deficiency. *New Engl. J. Med.* **2007**, *357*, 266–281. [[CrossRef](#)] [[PubMed](#)]
5. White, J.H. Regulation of intracrine production of 1,25-dihydroxyvitamin D and its role in innate immune defense against infection. *Arch. Biochem. Biophys.* **2012**, *523*, 58–63. [[CrossRef](#)] [[PubMed](#)]
6. Landrier, J.F.; Karkeni, E.; Marcotorchino, J.; Bonnet, L.; Tourniaire, F. Vitamin D modulates adipose tissue biology: possible consequences for obesity? *Proc. Nutr. Soc.* **2016**, *75*, 38–46. [[CrossRef](#)] [[PubMed](#)]
7. Landrier, J.F.; Marcotorchino, J.; Tourniaire, F. Lipophilic micronutrients and adipose tissue biology. *Nutrients* **2012**, *4*, 1622–1649. [[CrossRef](#)] [[PubMed](#)]
8. Bonnet, L.; Karkeni, E.; Couturier, C.; Astier, J.; Dalifard, J.; Defoort, C.; Svilar, L.; Martin, J.C.; Tourniaire, F.; Landrier, J.F. Gene expression pattern in response to cholecalciferol supplementation highlights cubilin as a major protein of 25(OH)D uptake in adipocytes and male mice white adipose tissue. *Endocrinology* **2018**, *159*, 957–966. [[CrossRef](#)]
9. Pereira-Santos, M.; Costa, P.R.; Assis, A.M.; Santos, C.A.; Santos, D.B. Obesity and vitamin D deficiency: A systematic review and meta-analysis. *Obes. Rev.* **2015**, *16*, 341–349. [[CrossRef](#)]
10. Vimalaswaran, K.S.; Berry, D.J.; Lu, C.; Tikkanen, E.; Pilz, S.; Hiraki, L.T.; Cooper, J.D.; Dastani, Z.; Li, R.; Houston, D.K.; et al. Causal relationship between obesity and vitamin D status: bi-directional Mendelian randomization analysis of multiple cohorts. *PLoS Med.* **2013**, *10*, e1001383. [[CrossRef](#)]
11. Bonnet, L.H.A.; Karkeni, E.; Couturier, C.; Astier, J.; Defoort, C.; Svilar, L.; Martin, J.C.; Tourniaire, F.; Landrier, J.F. Diet induced obesity modifies vitamin D metabolism and adipose tissue storage in mice. *J. Steroid. Biochem. Mol. Biol.* **2019**, *185*, 39–46. [[CrossRef](#)] [[PubMed](#)]

12. El-Khoury, J.M.; Reineks, E.Z.; Wang, S. Progress of liquid chromatography-mass spectrometry in measurement of vitamin D metabolites and analogues. *Clin. Biochem.* **2011**, *44*, 66–76. [[CrossRef](#)] [[PubMed](#)]
13. Grebe, S.K.; Singh, R.J. LC-MS/MS in the clinical laboratory—Where to from here? *Clin. Biochem. Rev.* **2011**, *32*, 5–31. [[PubMed](#)]
14. Mena-Bravo, A.; Ferreiro-Vera, C.; Priego-Capote, F.; Maestro, M.A.; Mourino, A.; Quesada-Gomez, J.M.; Luque de Castro, M.D. Quantitative analytical method to evaluate the metabolism of vitamin D. *Clin. Chim. Acta* **2015**, *442*, 6–12. [[CrossRef](#)] [[PubMed](#)]
15. Shah, I.; Petroczi, A.; Naughton, D.P. Method for simultaneous analysis of eight analogues of vitamin D using liquid chromatography tandem mass spectrometry. *Chem. Cent. J.* **2012**, *6*, 112. [[CrossRef](#)] [[PubMed](#)]
16. Fang, H.; Yu, S.; Cheng, Q.; Cheng, X.; Han, J.; Qin, X.; Xia, L.; Jiang, X.; Qiu, L. Determination of 1,25-dihydroxyvitamin D₂ and 1,25-dihydroxyvitamin D₃ in human serum using liquid chromatography with tandem mass spectrometry. *J. Chromatogr. B Anal. Technol. Biomed. Life Sci.* **2016**, *1027*, 1–26. [[CrossRef](#)] [[PubMed](#)]
17. Hedman, C.J.; Wiebe, D.A.; Dey, S.; Plath, J.; Kemnitz, J.W.; Ziegler, T.E. Development of a sensitive LC/MS/MS method for vitamin D metabolites: 1,25 dihydroxyvitamin D₂&3 measurement using a novel derivatization agent. *J. Chromatogr. B Anal. Technol. Biomed. Life Sci.* **2014**, 62–67. [[CrossRef](#)]
18. Abu Kassim, N.S.; Shaw, P.N.; Hewavitharana, A.K. Simultaneous determination of 12 vitamin D compounds in human serum using online sample preparation and liquid chromatography-tandem mass spectrometry. *J. Chromatogr. A* **2018**, *1533*, 57–65. [[CrossRef](#)]
19. Jenkinson, C.; Taylor, A.E.; Hassan-Smith, Z.K.; Adams, J.S.; Stewart, P.M.; Hewison, M.; Keevil, B.G. High throughput LC-MS/MS method for the simultaneous analysis of multiple vitamin D analytes in serum. *J. Chromatogr. B Anal. Technol. Biomed. Life Sci.* **2016**, *1014*, 56–63. [[CrossRef](#)]
20. Wan, D.; Yang, J.; Barnych, B.; Hwang, S.H.; Lee, K.S.; Cui, Y.; Niu, J.; Watsky, M.A.; Hammock, B.D. A new sensitive LC/MS/MS analysis of vitamin D metabolites using a click derivatization reagent, 2-nitrosopyridine. *J. Lipid Res.* **2017**, *58*, 798–808. [[CrossRef](#)]
21. Fabregat-Cabello, N.; Farre-Segura, J.; Huyghebaert, L.; Peeters, S.; Le Goff, C.; Souberbielle, J.C.; Cavalier, E. A fast and simple method for simultaneous measurements of 25(OH)D, 24,25(OH)₂D and the Vitamin D Metabolite Ratio (VMR) in serum samples by LC-MS/MS. *Clin. Chim. Acta* **2017**, *473*, 116–123. [[CrossRef](#)] [[PubMed](#)]
22. Xue, Y.; He, X.; Li, H.D.; Deng, Y.; Yan, M.; Cai, H.L.; Tang, M.M.; Dang, R.L.; Jiang, P. Simultaneous quantification of 25-hydroxyvitamin D₃ and 24,25-dihydroxyvitamin D₃ in rats shows strong correlations between serum and brain tissue levels. *Int. J. Endocrinol.* **2015**, 296531. [[CrossRef](#)]
23. Piccolo, B.D.; Dolnikowski, G.; Seyoum, E.; Thomas, A.P.; Gertz, E.R.; Souza, E.C.; Woodhouse, L.R.; Newman, J.W.; Keim, N.L.; Adams, S.H.; et al. Association between subcutaneous white adipose tissue and serum 25-hydroxyvitamin D in overweight and obese adults. *Nutrients* **2013**, *5*, 3352–3366. [[CrossRef](#)] [[PubMed](#)]
24. Lipkie, T.E.; Janasch, A.; Cooper, B.R.; Hohman, E.E.; Weaver, C.M.; Ferruzzi, M.G. Quantification of vitamin D and 25-hydroxyvitamin D in soft tissues by liquid chromatography-tandem mass spectrometry. *J. Chromatogr. B Anal. Technol. Biomed. Life Sci.* **2013**, *932*, 6–11. [[CrossRef](#)] [[PubMed](#)]
25. Ahonen, L.; Maire, F.B.; Savolainen, M.; Kopra, J.; Vreeken, R.J.; Hankemeier, T.; Myohanen, T.; Kylli, P.; Kostianen, R. Analysis of oxysterols and vitamin D metabolites in mouse brain and cell line samples by ultra-high-performance liquid chromatography-atmospheric pressure photoionization-mass spectrometry. *J. Chromatogr. A* **2014**, *1364*, 214–222. [[CrossRef](#)] [[PubMed](#)]
26. Burild, A.; Frandsen, H.L.; Poulsen, M.; Jakobsen, J. Quantification of physiological levels of vitamin D(3) and 25-hydroxyvitamin D(3) in porcine fat and liver in subgram sample sizes. *J. Sep. Sci.* **2014**, *37*, 2659–2663. [[CrossRef](#)] [[PubMed](#)]
27. Wang, Z.; Senn, T.; Kalhorn, T.; Zheng, X.E.; Zheng, S.; Davis, C.L.; Hebert, M.F.; Lin, Y.S.; Thummel, K.E. Simultaneous measurement of plasma vitamin D(3) metabolites, including 4beta,25-dihydroxyvitamin D(3), using liquid chromatography-tandem mass spectrometry. *Anal. Biochem.* **2011**, *418*, 126–133. [[CrossRef](#)] [[PubMed](#)]
28. Shah, V.P.; Midha, K.K.; Findlay, J.W.; Hill, H.M.; Hulse, J.D.; McGilveray, I.J.; McKay, G.; Miller, K.J.; Patnak, R.N.; Powell, M.L.; et al. Bioanalytical method validation—A revisit with a decade of progress. *Pharm. Res.* **2000**, *17*, 1551–1557. [[CrossRef](#)]

29. Armbruster, D.A.; Pry, T. Limit of blank, limit of detection and limit of quantitation. *Clin. Biochem. Rev.* **2008**, *29*, S49–S52.
30. Carrelli, A.; Bucovsky, M.; Horst, R.; Cremers, S.; Zhang, C.; Bessler, M.; Schrope, B.; Evanko, J.; Blanco, J.; Silverberg, S.; et al. Vitamin D storage in adipose tissue of obese and normal weight women. *J. Bone Miner. Res.* **2017**, *32*, 237–242. [[CrossRef](#)]



© 2019 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (<http://creativecommons.org/licenses/by/4.0/>).